

PLUG ELITE-SPK® ON THE INSERTER

- ▶ Insert a temporary shuttle suture as Vicryl suture size or a PDS USP 0 inside the rear eyelet of the device Elite-SPK®.


TIP! The shuttle sutures don't contribute to the repair stability and are removed once the construct is completed.

- ▶ Advance the mobile hook acting on the knob of the instrument. The lateral pin is out.

- ▶ Plug the the rear eyelet of the device Elite-SPK® through the mobile hook of the inserter.

⚠ WARNING! Use caution and insert the device in the correct direction: use the laser marking as reference. The device and inserter cannula must be aligned.

- ▶ Withdraw the knob, lock the system with the lateral pin and load the sutures in the distal eyelet of Elite-SPK®.

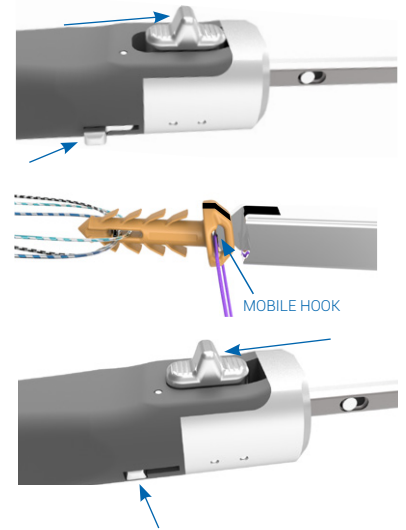

PERFORM THE REPAIR WITH ELITE-SPK®

- ▶ Tie each suture limb to the shuttle suture passing through the transosseous tunnel.

TIP! In presence of multiple shuttles, pull the medial limbs one by one for a easier transit through the transosseous tunnel.


It is recommended to use a para-acromial portal to drag the shuttle suture and the high-strength sutures more easily through the transosseous tunnel.

- ▶ Hold Elite-SPK® coaxial to the lateral access of the transosseous tunnel and draw the sutures to gently carry the implant forward through the deltoid.

TIP! Apply slight pressure or tension at the rear of the inserter to help the passage of the implant.

- ▶ Place Elite-SPK® on the lateral cortex. If it is done, use the slot previously performed with the punch to guide the input.
- ▶ Insert completely the body of Elite-SPK® inside the humerus and tighten the supporting underhead of the device to the lateral cortex.


TIP! A gently hammering might be required.

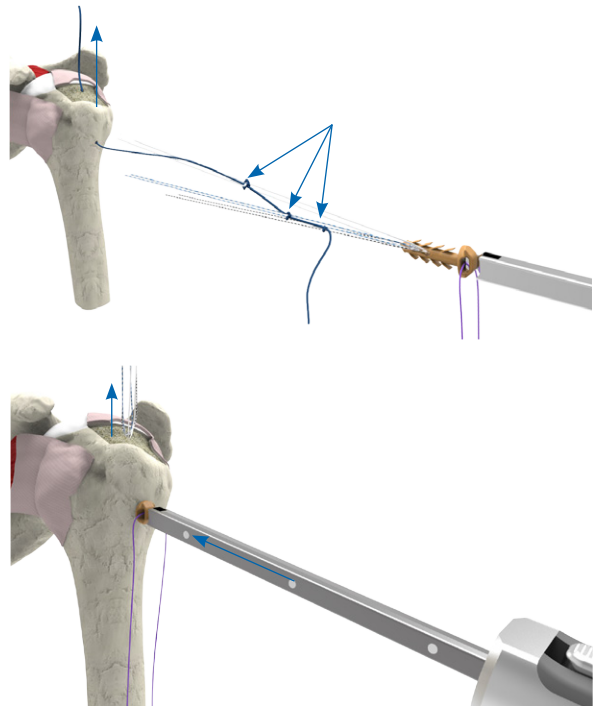

⚠ WARNING! Take care not to sink the implant into the cortex!

- ▶ Unlock the lateral pin of the inserter and pull the knob forward to release Elite-SPK®, then slide the inserter laterally to remove it.

- ▶ Suture the rotator cuff in the manner deemed most appropriate.

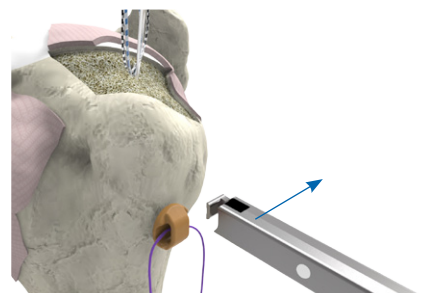

PLUG SHARC-FT® ON THE INSERTER

- ▶ Insert a temporary shuttle suture as Vicryl suture size or a PDS USP 0 inside the rear eyelet of the device Sharc-FT®.


TIP! The shuttle sutures don't contribute to the repair stability and are removed once the construct is completed.

- ▶ Advance the mobile hook acting on the knob of the instrument.

- ▶ Plug the the rear eyelet of the device Sharc-FT® through the mobile hook of the inserter.

WARNING! Use caution and insert the device in the correct direction: use the laser marking as reference. The device and inserter cannula must be aligned.

- ▶ Release the knob to set the implant on the inserter and load the sutures in the distal eyelet of Sharc-FT®.

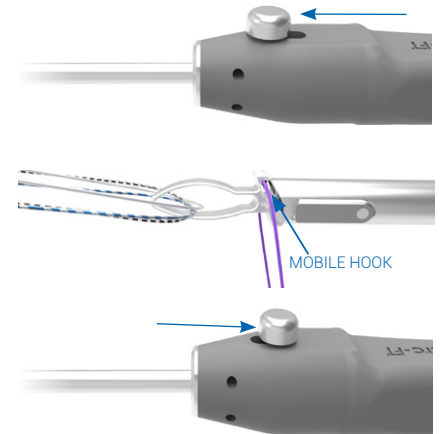

PERFORM THE REPAIR WITH SHARC-FT®

- ▶ Tie each suture limb to the shuttle suture passing through the transosseous tunnel.

TIP! In presence of multiple shuttles, pull the medial limbs one by one for a easier transit through the transosseous tunnel.

It is recommended to use a para-acromial portal to drag the shuttle suture and the high-strength sutures more easily through the transosseous tunnel.


- ▶ Hold Sharc-FT® coaxial to the lateral access of the transosseous tunnel and draw the sutures to gently carry the implant forward through the deltoid.

TIP! Apply slight pressure or tension at the rear of the inserter to help the passage of the implant.

- ▶ Place Sharc-FT® on the lateral cortex. If it is done, use the slot previously performed with the punch to guide the input.

- ▶ Insert completely the body of Sharc-FT® inside the humerus and tighten the supporting underhead of the device to the lateral cortex.

TIP! A gently hammering might be required.

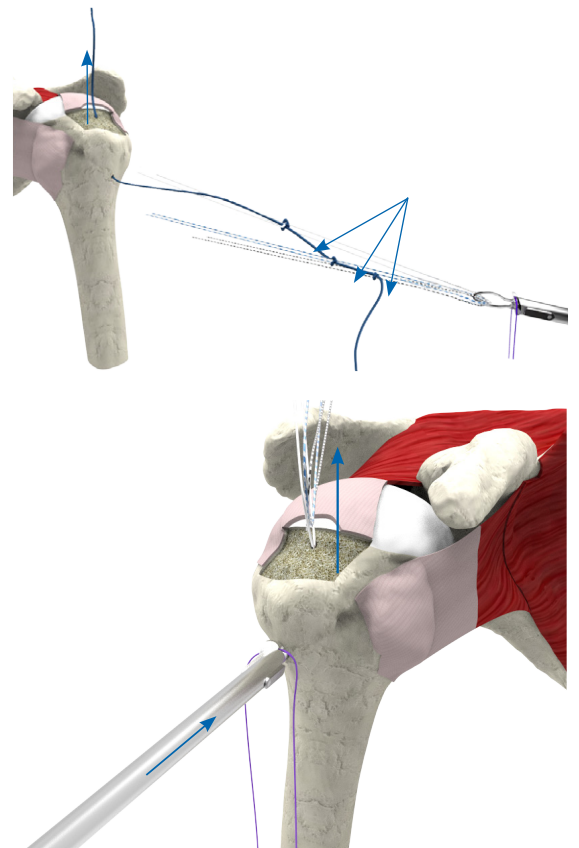

⚠ WARNING! Take care not to sink the implant into the cortex!

- ▶ Pull the knob forward to release Sharc-FT®, then slide the inserter laterally to remove it.

- ▶ Suture the rotator cuff in the manner deemed most appropriate.

